




The new modern style of the ERGOLINE offers a high technology look as well as with a hotel design providing all equipment for care and comfort.

Thick aluminium profile to resist to frontal shocks (patient bed or other heavy mobile equipment). Possibility to integrate a Multimedia movable arm for communication.

Also available for 2 beds on the same control trunking allowing total access for the room.

PROFIL

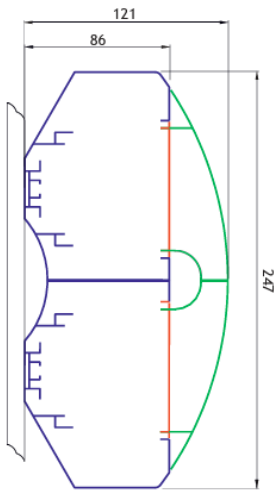

Section 247 mm x 121 mm

TOTAL MODULARITY



- Modular design of all components. For eventual modification on site clipping and unclipping cover system without cutting the housing.
- Possible flush mounting system (hygiene).
- Special insulated plastic housing for gas and electricity.

MEDICAL GAS OUTLETS



All types of gases outlets can be integrated with robustness and strongness (fixation in stainless steel bracket with aluminium cover painted white or grey).

Manufacturer medical gases outlet cover can be integrated as well as Biolume one.

The ventilation of compartments is assured by invisible baffle under the covers, for a better asepsis (in conformity with EN ISO 11197 Standards).

ELECTRICAL EQUIPMENT

- Bacteria resistant electrical sockets. Total adaptation to all international standards.
- Easy clip terminals wiring connections.

LIGHTINGS


- **General lighting** by 2x24 w fluocompact lamp. Polycarbonate smooth robust diffuser of 12 mm thickness.
- **Reading light** by 1x11w fluocompact lamp on rubber sheathed flexible arm.
- **Care light** by an extra 1x24 w integration into the general lighting controlled by switch on the unit.
- **Night light** : LED light integrated at the bottom of the panel. Discreet, efficient and easy maintenance.



Ambient lighting




Reading light

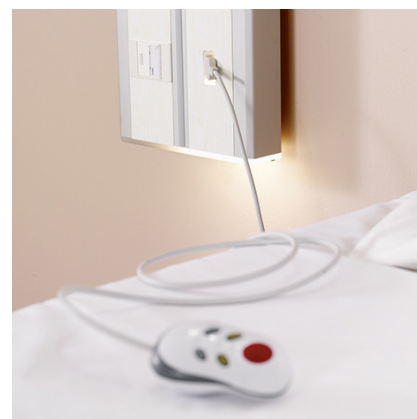

Night light

Ergoline									
	A	B	C	D	E	F	G	H	I
Ambiance 2x24w PL	145	135	105	80	85	75	170	165	140
Eclairage de soins* 2x24w PL + 1x24w PL + lampe de lecture 11w	540	305	190						

* sans la lampe articulée

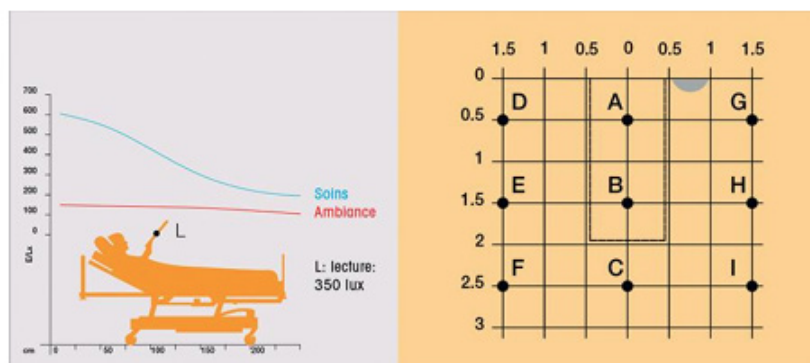

Eclairage direct	1 + 2x24w PL	Equipements électroniques HF
Lampe articulée	1 x 13w PL	Equipements électroniques HF
Eclairage veille	LED 1w	Convertisseur électronique HF

FURNITURE - SLIDING COVER



- Sliding cover to hide medical gases outlets and electrical sockets, in order to forget «hospital environment».
- The panel can be equipped with various furniture : Shelf, drawer, book support medical rail, infusion tube support, etc...
- This ergonomic configuration allows more space for care considering that the bedside is close to the bed (loading capacity of the shelf is about 20 kg).

TOP FINISH

- Peinture poudre époxy de couleur blanche, RAL 9016, finition mate et gris aluminium.
- Autres teintes possibles (sur demande).
- Embout fixé par clipage intérieur
- Couvertures décor bois (option), nuancier sur demande

TECHNICAL CHARACTERISTICS

- Lighting systems in conformity with Standard N F EN 12464-1 and French Lighting association (recommendation, Feb 2000)
- Electrical registration Class I NF EN 60601-1.
- IP20 Rating protection.
- 850°C incandescent wire rating (NF EN 60695-2-11).
- Cardboard and shrink wrap packaging.
- Epoxy paint finish white RAL 9016 (or grey RAL 9006).
- Materials : Aluminium and stainless steel and polycarbonate.
- Average weight: 40 kg/ lm.
- Suitable wall mounting plugs and screws.
- Easy cleaning and disinfection by soft cloth.

